## PROPLUS<sup>®</sup> 2650KGWA

**ISO13485:2003  
CERTIFIED**

**USER INSTRUCTIONS  
INSTRUCCIONES DE USO  
MODE D'EMPLOI  
INSTRUÇÕES DE USO  
BEDIENUNGSANLEITUNG  
INSTRUCTIES VOOR DE GEBRUIKER  
BRUGSVEJLEDNING  
KÄYTTÖOHJEET  
ISTRUZIONI PER L'USO  
BRUKSANVISNING  
INSTRUKCJA DLA UŻYTKOWNIKA  
BRUKSANVISNING  
用户说明**

**PELSTAR, LLC 11800 Austin Ave. Alsip, IL 60803**

**Distributed by:**

**WelchAllyn<sup>®</sup>**

**[www.welchallyn.com](http://www.welchallyn.com)**

Patents 7550682B2 and D508655 (US); Patents 0149984/1-8 and Pending (EU);  
Patents ZL200480031711X and 200430004551.2 (PRC) © Pelstar, LLC 2010

Rev. 20101108

---

## **PROPLUS<sup>®</sup> 2650KGWA**

### **TABLE OF CONTENTS**

<b>ENGLISH .....</b>	<b>2</b>
<b>ESPAÑOL .....</b>	<b>14</b>
<b>FRANÇAIS-.....</b>	<b>26</b>
<b>PORTUGUÊS .....</b>	<b>38</b>
<b>DANSK .....</b>	<b>50</b>
<b>DEUTSCH.....</b>	<b>62</b>
<b>NEDERLANDS .....</b>	<b>74</b>
<b>SUOMI .....</b>	<b>86</b>
<b>ITALIANO .....</b>	<b>98</b>
<b>NORSK .....</b>	<b>110</b>
<b>POLSKI .....</b>	<b>122</b>
<b>SVENKA .....</b>	<b>134</b>
<b>简体中文 .....</b>	<b>146</b>

**PROPLUS<sup>®</sup> 2650KGWA**

**DIGITAL WHEELCHAIR PLATFORM SCALE  
USER INSTRUCTIONS**

**Thank you** for your purchase of this Health o meter<sup>®</sup> Professional product.  
Please read this manual carefully and keep it for easy reference or training.

**TABLE OF CONTENTS**

CAUTIONS AND WARNINGS ..... 3

SPECIFICATIONS..... 3

CERTIFICATIONS / CONNECTIVITY / DISPOSAL / SERVICE ..... 4

ASSEMBLY INSTRUCTIONS..... 5

SET UP / BATTERIES..... 6

QUICK START INSTRUCTIONS ..... 7

MAINTENANCE / CALIBRATION ..... 9

TROUBLESHOOTING..... 10

DETAILED VIEW OF SCALE ..... 11

PARTS FOR MAINTENANCE ..... 12

WARRANTY ..... 13

## CAUTIONS AND WARNINGS

**To prevent injury and damage to your scale, please follow these instructions very carefully.**

- Do not transport the scale while the platform is loaded. Unplug the adapter and store both the adapter and cord before moving the scale. Be sure the scale is tilted enough to assure it is rolling properly on its wheels when it is moved for cleaning or transported from one area to another.
- For accurate weighing the scale must be placed on a flat, stable surface.
- For accurate weighing, verify proper operation according to the procedure described in this manual before each use. Do not exceed recommended weight limit of 400kg for this scale.
- Do not use in the presence of flammable materials.
- Operating at other voltages and frequencies than specified could damage the equipment.
- If the "LOW BAT" indicator activates, for accurate weighing, replace the batteries or connect the scale to a power source, using the power adapter, as soon as possible, depending on use—see Connectivity section.

## SPECIFICATIONS

### GENERAL

The ProPlus® 2650KGWA Wheelchair Platform Digital Scale for Patients uses highly sophisticated microprocessor technology. Each precision scale is designed to provide accurate, reliable and repeatable weight measurements and features that make the weighing process simple, fast and convenient.

These scales are set up to use motion-sensing weighing technology, to determine the actual weight of a moving patient. The scales may be changed to measure live weight (see the ProPlus® Display Module User Manual for instructions on how to change the scale settings).

Weight is displayed in kilograms. The scale features a remote display which may be held by the operator or placed on a table, desk, wall or other flat surface for privacy and convenience in use. The 2650KGWA has 2 wheels built in for easy mobility. The scale can be operated using a power adapter (included) or by 6-D cell batteries (not included), depending on use—see Connectivity section, below.

### 2650 KGWA SPECIFICATIONS

<b>Capacity and Resolution</b>	<b>2650KGWA,:</b> Max 400kg / Min 1kg, e=0.1kg (Class III)	
<b>Power Requirements</b>	120V-240V AC, 50-60Hz - 9V DC adapter, ADAPTER MODEL ADPT30 (INCLUDED) or 6 D-cell batteries (not included)	
<b>Environmental</b>	Operating temperatures: 10°C to 35°C (50°F to 95°F) Storage temperatures: 0°C to 50°C (30°F to 125°F) Maximum Humidity: 85% RH	
<b>Physical Dimensions</b>	<b>Platform - 2650 Series</b> Length: 82.6 cm (32-½") Width: 91.4 cm (36") Height: 6.2 cm (2-½")	<b>Assembled Size and Weight</b> Length: 105.4 cm (41-½") Width: 143 cm (56-¼") Height: 8 cm (3-⅛") Weight: 36.5 kg (80.5 lb)

# CERTIFICATIONS / CONNECTIVITY / DISPOSAL

## CERTIFICATION DESCRIPTIONS

	<p>Scales carrying this symbol meet the Precision Class III standards of the Directive 90/384 EEC</p>
	<p>Products carrying this symbol meet the requirements of:</p> <ol style="list-style-type: none"> <li>1. Directive for non-automatic weighing instruments 90/384/EEC</li> <li>2. Directive governing medical devices 93/42/EEC</li> </ol>
	<p>Pelstar, LLC has been officially certified as a manufacture of Health o meter® Professional medical devices to be in compliance with stringent legal requirements. The Pelstar, LLC quality assurance system covers the development, production, sales and service of medical scales and measuring systems.</p>

## 2650KGWA SCALE CONNECTIVITY INFORMATION

These scales are designed to connect to a computer, monitor or other devices via USB or serial ports, to allow reliable transmission of weight, height and BMI from the scale. This connectivity helps close a source of potential medical errors caused by incorrectly copying and then recording weight information taken. This scale is EMR-EHR compatible. For more information regarding how to connect this scale to devices please see the ProPlus® Display Module User Manual for instructions or contact Customer Service at one of the numbers below.

## 2650KGWA SCALE DISPOSAL

This Health o meter® Professional scale must be disposed of properly as electronic waste. Follow the national, regional or local regulations which apply to you for disposal of electronic waste or batteries. Do not dispose of this device in domestic waste.

## CUSTOMER SERVICE INFORMATION

For further information or telephone support, please contact Welch Allyn Customer Service at one of the numbers listed below

- European Call Center - (+353) 46-90-67790
- France - (+33) 1-55-69-58-49
- Germany - (+49) 695-098-5132
- Italy - (+39) 026-968-2425
- Netherlands- (+31) 202-061-360
- Spain - (+34)917-499-357
- Sweden - (+46) 85-853-6551
- United Kingdom - (+44) 207-365-6780
- Latin America - (+1) 305-669-9003
- Middle East - (+49) 7477-927186
- South Africa - (+27) 11-777-7555
- Australia - (+61) 2-9638-3000
- China - (+86) 21-6327-9631

# 2650KGWA ASSEMBLY INSTRUCTIONS

## **Before Assembly (2650KGWA):**

The ProPlus® 2650KGWA Wheelchair Patient Digital Platform Scale is shipped in one carton. Carefully inspect the carton for shipping damage before unpacking. If damage is found, contact Customer Service. The following information details what you will find inside the main carton.

**Remove each** assembly from the carton and unwrap packing material carefully, to prevent scratching the unit's parts. Set the carton aside, for storage. To avoid damaging the scale parts when unpacking it, **do not use a box cutter, knife, scissors or any sharp object to open the protective inner packaging.**


## **Parts List:**

- (1) Base / Platform Assembly (with 2 ramps, handle, wheels attached)
- (1) User information
- (1) Box Containing:
  - (1) Display Assembly with attached Load Cell Connector Cable
  - (1) Small Box Containing:
    - (1) Power Adapter
    - (4) Feet

## **Assembling a 2650KGWA Scale:**

**This scale takes only minor assembly.**

1. Place scale base on a flat, level, dry floor—place upside down carefully.
2. Thread the 4 feet into the underside of the base. Each foot should be completely hand tightened and then backed out 1-½ full turns. This allows for proper load cell operation for accuracy of the scale.
3. Carefully tilt the scale back to an upright position on the floor with care to watch the ramp—it may flip down when the scale is tipped back. **DO NOT DROP THE SCALE.** Adjust the feet, to level the scale. Be certain that the feet are **NOT** completely tightened into the base (see Step 2, above).
4. Open the scale to operational position.
5. Connect the load cell cable connector to the load cell port on the scale base.
6. Connect the power adapter cord to the power port in the display, if the adapter will be used. The cord will be outside of the display, as shown in Figure 1. It is optional to run the cord through the channel in the back of the display, if you prefer. To accomplish this, remove the 4 screws (A) in the back with a Phillips screwdriver, remove the cover, run the power cord (B) through the channel, and then replace the cover, secure it with the 4 screws and tighten.
7. The scale is now assembled. Plug the scale adapter into a power source and / or insert batteries—see the batteries section / Figure 2 for more information regarding batteries.

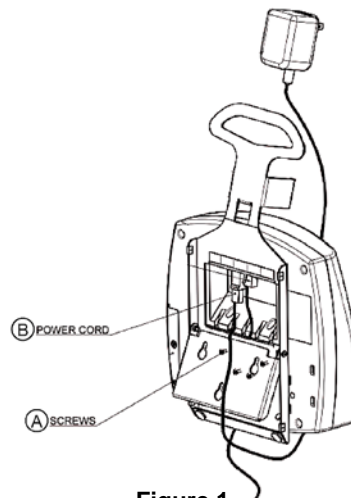**Figure 1**

## SET UP / BATTERIES

ProPlus® scales are shipped with the Menu Lock function enabled. The Menu Lock function disables certain buttons on the scale's keypad, preventing unauthorized users from making changes to your scale settings. The Menu Lock can be bypassed temporarily to allow you to customize your scale with the functionality you prefer. Please refer to the ProPlus® Display Module User Manual for instructions on how to bypass the Menu Lock.

### Preparing the Scale for Use


1. Remove protective plastic films covering the keypad and display.
2. Place batteries in the battery holder (see "Batteries" information in Figure 2, below).
3. Plug the scale's power adapter into the power source.
4. Be sure all objects are removed from the scale's platform, and then press the **ON/OFF** or **I/O** button to turn the scale on. The display will show "Health o meter ProPlus" and then "0.0KG".

**Note:** To adjust the display backlight and/or contrast, please refer to the ProPlus® Display Module User Manual.

**Note:** The scale will always default to the settings and units last used.

5. When you disconnect the scale's power adapter from the power source, the scale will switch to battery power, if batteries are installed. **Note:** If the set up procedure failed, refer to the troubleshooting instructions. If the problem is not corrected, refer to qualified service personnel.

### Replacing / Inserting Batteries

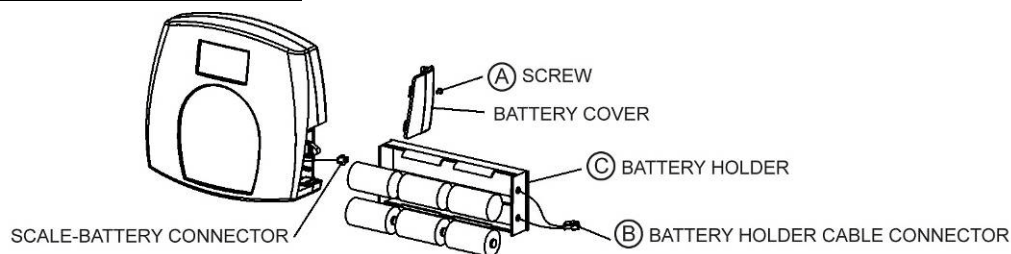**Figure 2 - Batteries**

1. Unplug the scale.
2. Remove the battery cover from the display assembly—refer to Figure 2 (A) above.
3. Disconnect battery holder cable connector from the scale-battery connector (B).
4. Carefully remove the battery holder (C) by sliding it out of the display assembly.
5. Replace the batteries with new ones. We recommend EVEREADY Energizer® e2™ batteries.
6. Carefully slide the battery holder (C) into the display assembly.
7. Plug the battery holder cable connector to the scale-battery connector (B).
8. Attach the battery cover to the display assembly and install the screw (A).

**Warning:** If the scale will not be used for some time, remove batteries to avoid a safety hazard. Disposal of batteries should be done according to the national, regional or local regulations which apply to you for disposal of batteries.

## SET UP / TRANSPORT

### Preparing the Scale for Transport

The scale can easily be folded for movement from one room to another and includes pins on one ramp, to carry the indicator. To prepare the scale for transport, follow these steps:

1. Fold both ramps--left and right (A) onto the platform while the scale is on the floor (see Figure 3).
2. Place the indicator (B) onto the pins on one of the folded ramps by sliding the slots that are on the back of the indicator onto the raised pins.

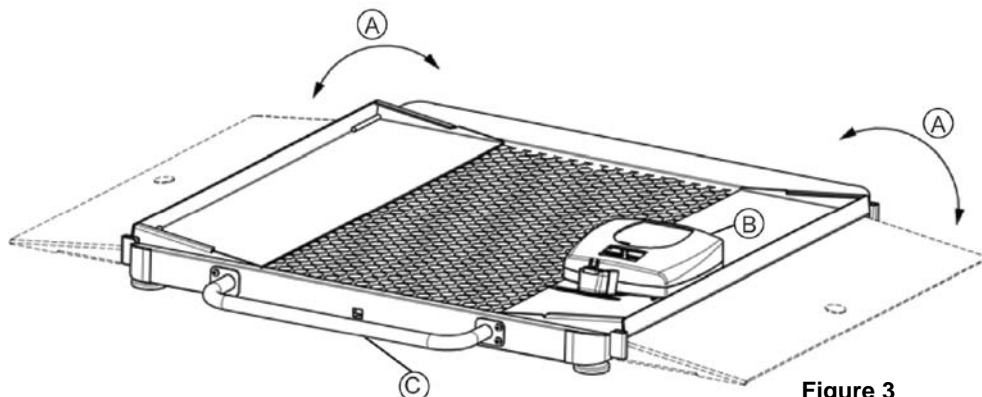
Figure 3

3. Lift the scale onto its side by pulling up on the handle (C). The scale will now be ready to transport, rolling on its wheels (see Figure 4).

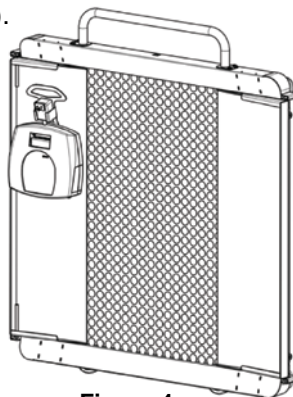
Figure 4

## QUICK START INSTRUCTIONS

**Additional, complete instructions regarding all of the display and keypad functions and menus can be found in the ProPlus® Display Module User Manual included with this scale.**

**NOTE:** This scale will always default to the settings and units last used, when it is turned on.

### WEIGHING A PATIENT

1. Make sure there is no object on the weighing platform.
2. Press the **ON/OFF** or **I/O** button to turn the scale on.
3. Wait until "0.0Kg" and "ZERO" on the left side of the display appear.
4. Roll the patient onto the scale. The scale will calculate the patient's weight but the display will not show a weight until the stable weight of the person is determined. Depending on the movement of the patient on the scale, it may take several seconds for the scale to lock onto the weight.
5. Once a weight has been displayed, while the patient is still on the scale, we recommend performing a second weighing by pressing the REWEIGH button.
6. After the weight measurement has been taken and is recorded, roll the patient roll off the scale.
7. Press the **ON/OFF** or **I/O** button to turn the scale off.

## QUICK START INSTRUCTIONS (CONTINUED)

### CALCULATING BMI

Body Mass Index is another name for Quetelet Index, a formula that is a convenient and reliable indicator of obesity. **NOTE:** The scale will not calculate BMI for a patient that is 12 kg or less. If a height is not entered for the patient within 30 seconds of pressing the BMI button, the scale will return to the normal weighing mode.

1. Complete steps 1 to 4 above for "Weighing A Patient".
2. After obtaining the patient's weight on the scale, press the BMI (2) button while the scale is on.
3. The scale will prompt you to enter the patient's height. Use the keypad to enter the height.  
**NOTE:** When weighing in kilograms, the height must be entered in 0.1 centimeter increments.
4. Once height is keyed in, press **ENTER**.
5. The display will show the patient's **BMI** (Body Mass Index).
6. Press **EXIT** to return to the normal weighing mode of the scale.

### TARE FUNCTION

When using the scale, the weight of an object such as a wheelchair can be subtracted from the total weight to determine the weight of the patient alone. The Tare Function automatically performs this subtraction. ProPlus® scales allow tare weight to be entered manually by using the keypad or automatically, as outlined in this section.

#### AUTOMATIC (PUSH BUTTON TARE)

1. Wait until "0.0Kg" and "ZERO" appear on the left side of the display.
2. With all weight off the scale, press the **ON/OFF** or **I/O** button to turn on the scale.
3. Place the wheelchair on the scale. The display will show a value for the wheelchair's weight.
4. Press the TARE (9) button. The scale will zero out and the word "TARE" will appear in the upper right side of the display to indicate there is a tare value for this patient's wheelchair.
5. Remove the weight from the scale. The display will show a negative weight equal to the weight of the wheelchair, once it is removed from the scale.
6. Place the patient and wheelchair on the scale. The scale will automatically deduct the tared weight of the wheelchair from the gross weight of the patient and wheelchair together.
7. The weight of the patient will appear on the scale.
8. The tared value is stored in memory until changed, cleared or the scale is turned off.
9. To remove the Tare, press the TARE (9) button again. The word "TARE" will disappear from the display and the tare value will be deleted from memory. You can also press the EXIT button. The scale will resume normal operation.

#### MANUAL (KEYBOARD OR KEYPAD TARE)

1. Make sure there is no object on the weighing platform.
2. Press the **ON/OFF** or **I/O** button to turn on the scale.
3. Wait until "0.0Kg" and "ZERO" appear on the left side of the display.
4. Press the **TARE (9) button**. The user will be prompted to enter the tare value.
5. Use the keypad to enter the weight of the shoes or other object and then press **ENTER**.
6. The value entered will show as a negative value on the display.
7. Place the patient and the tared object on the scale. The scale will automatically deduct the tared weight, as entered, from the gross weight of the patient with the tared object.
8. The weight of the patient will appear on the scale.
9. The tared value is stored in memory until changed, cleared or the scale is turned off.
10. To remove the Tare, press the **TARE (9) button** again. The word "TARE" will disappear from the display and the tare value will be deleted from memory. You can also press the **EXIT** button. The scale will resume normal operation.

#### REMOVING THE TARE VALUE

Press the **TARE (9) button** again to delete a tare value from memory. The word "TARE" will disappear from the display and the tare value will be deleted from memory. You can also press the **EXIT** button. The scale will resume normal operation.

## MAINTENANCE

The following pages provide instructions for maintenance, cleaning, troubleshooting, and operator replaceable parts for the ProPlus® 2650KGWA scale. Maintenance operations other than those described in this manual should be performed by qualified service personnel.

### **MAINTENANCE**

Before first use and after periods of non-use, check the scale for proper operation and function. If the scale does not operate correctly, refer to qualified service personnel.

1. Check overall appearance of the total scale for any obvious damage, wear and tear.
2. Inspect the power adapter for cord cracking or fraying or for broken or bent prongs.

### **CLEANING**

Proper care and cleaning is essential to ensure a long life of accurate and effective operation. **Disconnect the scale from the adapter power source.**

1. Clean all external surfaces with a clean damp cloth or tissue. Mild antimicrobial or antibacterial soap and water solution may be used. Dry with a clean soft cloth.
2. Do not immerse the scale into cleaning or other liquid solution. Do not gas or steam sterilize the scale, nor place any part of it through an autoclave.
3. Do not use Isopropyl Alcohol or other solutions to clean the display surface.
4. Do not use abrasive cleaners.

## 2650KGWA CALIBRATION

**This scale has been factory calibrated and does not require calibration prior to use.** For maintenance, the scale can be calibrated. Calibration of your ProPlus® scale is performed using kilos (kg).


Operator Action	Display
1. Make sure the display is turned off, then press and hold the HOLD/RELEASE button while pressing the <b>ON/OFF</b> or <b>I/O</b> button.	“Enter load weight: “200.0KG”
2. Using the ▲▼ keys and/or the numeric keypad, set the calibration load to at least half the capacity of the scale (for the best results and greater accuracy, use the full capacity of the scale.	“Enter load weight: “200.0KG”
3. Make sure there is nothing on the scale, then press the ENTER button.	Zero calibration: Please clear the scale
4. Please wait 1-3 seconds, without touching or shaking the scale, to allow the Zero Calibration to finish.	Zero calibration: Please wait
5. No action required by operator.	Weight calibration: Put: XXX.X
6. Load the scale with the required weight and press the ENTER button. Please wait 1-3 seconds, without touching or shaking the scale, to allow the Calibration to finish.	Weight calibration Please wait
7. If calibration is successful, the display will show “Factor OK!”. If the calibration process failed, it will show “Factor Fail!”.	Factor Update... “Factor OK!” or “Factor Fail!”
8. If successful, the display will begin to reboot.	Rebooting... Please wait
9. The display will return to the start up message.	Health o meter ProPlus
10. Please wait until the scale resumes normal operation and then remove the weight from the scale.	Health o meter ProPlus

# TROUBLESHOOTING

Refer to the following instructions to check and correct any failure before contacting Customer Service.

<b>SYMPTOM</b>	<b>POSSIBLE CAUSE</b>	<b>CORRECTIVE ACTION</b>
Scale does not turn on	<ol style="list-style-type: none"> <li>1. Dead Battery</li> <li>2. Faulty electrical outlet</li> <li>3. Bad power supply</li> </ol>	<ol style="list-style-type: none"> <li>1. Replace batteries</li> <li>2. Use a different outlet</li> <li>3. Replace adapter</li> </ol>
Questionable weight or the scale does not zero	<ol style="list-style-type: none"> <li>1. External object interfering with the scale</li> <li>2. The display did not show "0.0" before weighing</li> <li>3. Scale is not placed on a level floor</li> <li>4. Scale is out of calibration</li> <li>5. Improper Tare</li> </ol>	<ol style="list-style-type: none"> <li>1. Remove interfering object from the scale</li> <li>2. Ask the patient to get off the scale, zero the scale and begin weighing process again</li> <li>3. Place the scale on a level floor and begin weighing process again</li> <li>4. Check weight with known weight value</li> <li>5. Place the item to be tared on the scale. Press REWEIGH. Once the weight of the item is displayed, press TARE. Place the patient and the tared item back on the scale. Press REWEIGH again</li> </ol>
Weighing is performed but the display shows "weigh" and "reweigh" every few seconds; the weighing process takes too long and no weight is displayed.	The patient is not standing still	Ask the patient to stand still or you can change to live weight setting
The display shows "Overload" message	The load on the scale exceeds the capacity (400 kg)	Remove the excess weight and use the scale according to its limits
The display shows "LOW BAT" message	The batteries are empty	Replace batteries according to instructions
The display shows "Load Cell Error" message	There is a problem with one or more load cells or the load cell cable is disconnected.	Check load cell cable connection at the display and platform assembly ports. If the problem is not corrected, refer to qualified service personnel to replace the defective load cell

DETAILED VIEW OF 2650KGWA SCALE
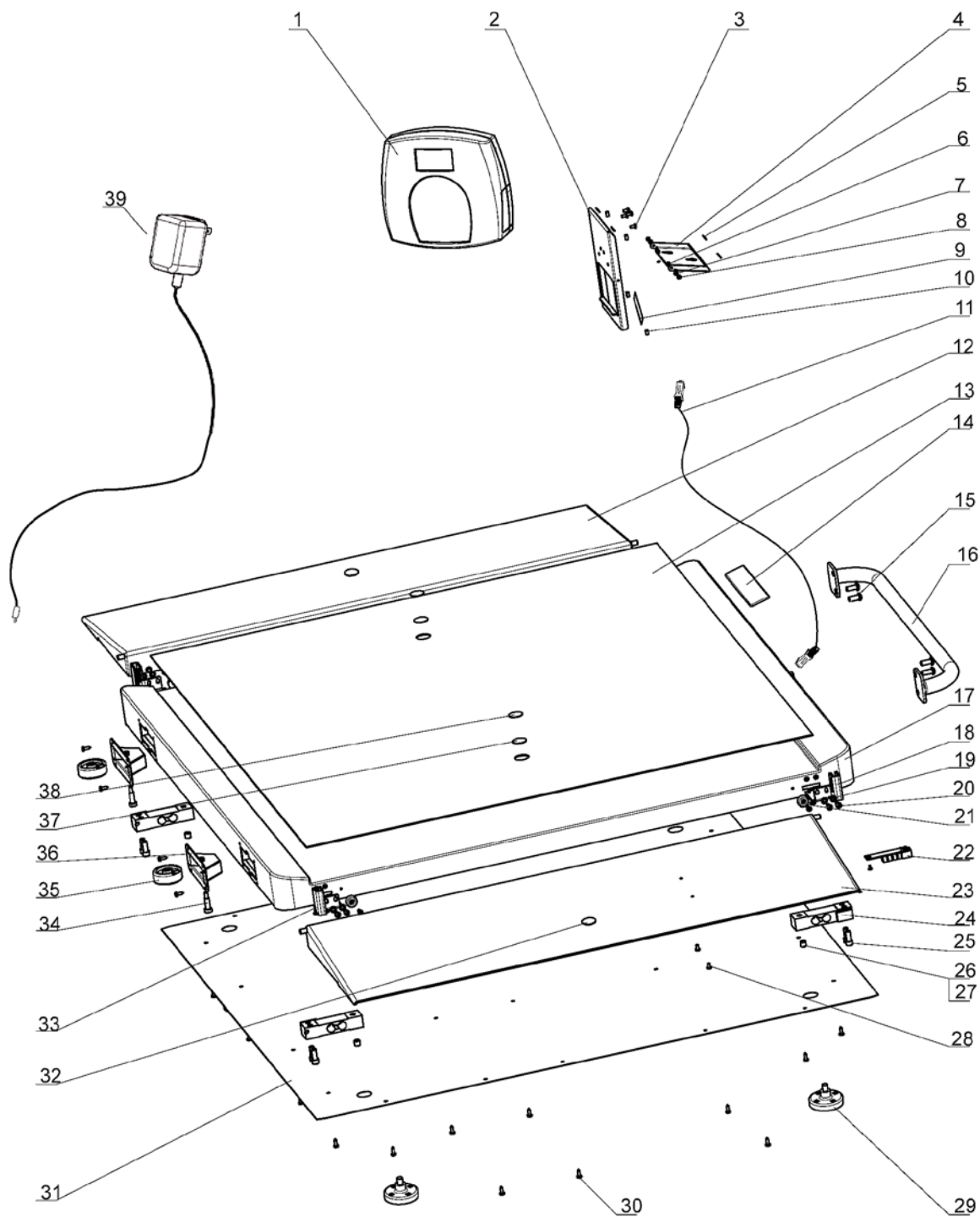
## 2650KGWA PARTS FOR MAINTENANCE

Field Replaceable Parts			
Key No.*	Part No.	Description	Qty.
1	B1109401-0	DISPLAY ASSEMBLY	1
11	B2033701-0	MAIN BOARD TO JUNCTION BOX LOAD CELL CABLE	1
18	B3823601-0	RIGHT RAMP BRACKET	2
22	B2138901-0	CONNECTOR BOARD	1
24	B1304202-0	LOAD CELL	4
29	B2266701-0	FOOT ASSEMBLY	4
33	B3823701-0	LEFT RAMP BRACKET	2
39	ADPT30	POWER ADAPTER 120V-240V	1
	PROKIT	INCLUDES ADAPTER, DISPLAY ASSY, 4 LOAD CELLS, CONNECTOR BOARD, 4 FOOT ASSEMBLY	1
Below listed items are for reference only (refer to detailed view on page 11)			
2		DISPLAY STAND	1
3		PHILLIPS PAN HEAD SCREW, M4	4
4		DISPLAY STAND FOOT	1
5		NON SKID PADS	4
6		SELF LOCK NUT, M3	2
7		WASHER, M3	2
8		SCREW, M3x10	2
9		INPUT / OUTPUT POSTS LABEL	1
10		BUMPER	4
12		RAMP #1	1
13		MAT	1
14		HEALTH O METER LABEL	1
15		SCREW, M8x19	4
16		HANDLE	1
17		PLATFORM ASSEMBLY	1
19		WASHER, M5	16
20		SCREW, M5x10	16
21		RAMP MAT	4
23		RAMP	1
25		SOCKET HEAD CAP SCREW 5/16-18x¾	8
26		SCREW, M10x10	4
27		LOAD CELL SAFETY SCREW	4
28		SCREW, M4x6	2
30		SHEET METAL SCREW, #8x½	26
31		BOTTOM COVER	1
32		SCOTCH MATE HOOK	2
34		SHOULDER BOLT	2
35		WHEEL	2
36		WHEEL CLIP ASSEMBLY	2
37		DOUBLE STICK TAPE	2
38		SCOTCH MATE LOOP	2

\* See page 11 for a detailed view of parts by number.

# WARRANTY

## 3- YEAR LIMITED WARRANTY

### What does the Warranty Cover?

Health o meter® Professional ProPlus® scales are warranted from date of purchase against defects of materials or in workmanship for a period of three (3) years. If you have an equipment problem that you cannot resolve, and the Customer Service provider determines that a defect of material or in workmanship exists, the customer's sole remedy will be repair or replacement of the scale at no charge. Replacement will be made with a new or remanufactured product or component. If the product is no longer available, replacement may be made with a similar product of equal or greater value. All parts including repaired and replaced parts are covered only for the original warranty period.

### Who is Covered?

The original purchaser of the product must have proof of purchase to receive warranty service. Please save your invoice or receipt. Authorized dealers selling Pelstar products do not have the right to alter, or modify or any way change the terms and conditions of this warranty.

### What is Excluded?

**Your warranty does not cover normal wear of parts or damage resulting from any of the following:** negligent use or misuse of the product, use on improper voltage or current, use contrary to the operating instructions, abuse including tampering, damage in transit, or unauthorized repair or alternations. Further, the warranty does not cover Acts of God, such as fire, flood, hurricanes and tornadoes. This warranty gives you specific legal rights, and you may also have other rights that vary from country to country, state to state, province to province or jurisdiction to jurisdiction.

**To get Warranty Service make sure you keep your sales receipt or document showing proof of purchase.** Call the Customer Service Center nearest you for assistance--please see page 4 of this manual. **Before returning a product, you must obtain authorization from Customer Service.**

Customer Service personnel will give you a Service Notification number. Please note this number on the outside of your shipping box. Returns without a Service Notification number will not be accepted for delivery. **Attach your proof of purchase to your defective product along with your name, address, daytime telephone number and description of the problem. Carefully package the product and send it with shipping and insurance prepaid to the address you are given by Customer Service.**

Write, "Attention SN# \_\_\_\_\_", and the shipping address, on the shipping label, following Customer Service instructions.

PRODUCT LIFE: 10 YEARS / 100,000 CYCLES

QNET BV  
HOMMERTERWEG 286  
6436 AM AMSTENRADE  
THE NETHERLANDS

PELSTAR, LLC

11800 South Austin Avenue - B • Alsip, IL 60803 • USA

www.homscales.com

ProPlus® is a registered trademark of Pelstar, LLC. Health o meter® is a registered trademark of Sunbeam Products, Inc. used under license. Health o meter® Professional products are manufactured, designed and owned by Pelstar, LLC and are distributed by Welch Allyn outside the US and Canada. We reserve the right to improve, enhance or modify Health o meter® Professional product features or specifications without notice.

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